

# FINAL REPORT

Assessment of Antimicrobial Activity

Using a Time-Kill Procedure

Order Number: 551713520

**PREPARED FOR**


Franke Kindred Canada | 1000 Franke Kindred Road  
Midland, Ontario L4R 4K9

Sneha Panchal, M.Sc., RMCCM

Microbiology Laboratory Manager

12/6/2017

EMSL Canada Inc.  
2756 Slough Street, Mississauga, ON L4T 1G3  
Phone: (289) 997-4602 Fax: (289) 997-4607 Web: [www.emsl.com](http://www.emsl.com)

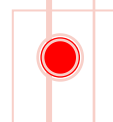

## CERTIFICATE OF ANALYSIS

---

**CLIENT:** FRANKE KINDRED CANADA

**PRODUCT:** OZONATED TAP WATER

**CONTACT:** CORY MACEY

**SAMPLE RECEIVED:** 11/16/2017

**PROJECT:**

**REPORT DATE:** 12/06/2017

ASSESSMENT OF ANTIMICROBIAL  
ACTIVITY USING ASTM METHOD  
E2315

**CHALLENGE BACTERIA:**

*Clostridium difficile*

### I. EXPERIMENTAL SUMMARY

The testing procedure was designed after discussions between EMSL Canada Inc. and Franke Kindred Canada. The procedure is based on ASTM E2315 method guidelines and conducted on two ozonated water samples of concentrations 2.0 ppm and 5.0 ppm, to demonstrate its effectiveness at killing *Clostridium difficile* (ATCC 9689). The testing was conducted in the Mississauga Microbiology Laboratory.

### II. PROCEDURE

The testing was done to determine the effectiveness of ozonated water samples of concentrations 2.0 ppm and 5.0 ppm (provided by Cory Macey) at killing *Clostridium difficile* (ATCC 9689) for 30 seconds and 15 minutes exposure times.

**Culture preparation:**

A pure culture of *Clostridium difficile* (ATCC 9689) was streaked on to Healthlink Anaerobic Blood Agar plates and incubated at 35°C for up to 72 hours under anaerobic conditions.

The control suspension of the test microorganism was standardized to a minimum concentration of  $1.0 \times 10^6$  CFU/mL, by dilution in a buffered saline solution.

Test and control substances were dispensed in identical volumes to sterile test tubes.

Independently, Test and Control substances were inoculated with the test microorganism, mixed and incubated.

Control suspensions were immediately plated to represent the concentration present at the start of the test, or time zero.

At the conclusion of each contact time, a volume of the liquid test solution was neutralized.



Dilutions of the neutralized test solution were plated on Healthlink Anaerobic Blood Agar plates and incubated at 35°C for up to 72 hours under anaerobic conditions to determine the surviving microorganisms at the respective contact times.

Reductions of microorganisms were calculated by comparing initial microbial concentrations to surviving microbial concentrations.

All tests were performed in duplicates and counts averaged.

**Calculations:**

Calculations were based on the following:

Percentage reduction =  $(B-A/B) \times 100$

Log 10 Reduction =  $\text{Log}(B/A)$

Where:

B = Number of viable test microorganisms in the control substance immediately after inoculation

A = Number of viable test microorganisms in the test substance after the contact time

**III. EXPERIMENTAL RESULTS**

Test Microorganism	Test Substance	Contact Time	CFUs/mL	Percent Reduction Compared to Control at Time Zero	Log10 Reduction compared to Control at Time Zero
<i>C.difficile</i>	CONTROL	Time Zero	3.6E+04	n/a	
	<b>SAMPLE A:</b> 2.0 PPM CONC.	30 seconds	2.8E+03	92.22%	1.11
Test Microorganism	Test Substance	Contact Time	CFUs/mL	Percent Reduction Compared to Control at Time Zero	Log10 Reduction compared to Control at Time Zero
	<b>SAMPLE B</b> 5.0 PPM CONC.	30 seconds	1.1E+03	96.94%	1.51
		15 mins	<1.00E+01	>99.98%	>4.90



#### IV. CONCLUSIONS/OBSERVATIONS

Sample B (concentration 5.0 ppm) caused a 99.98% reduction on C.difficile with 15 minutes exposure. Sample B also proved to be more effective for the 30 seconds exposure time compared to Sample A for the same time point.



---

Sneha Panchal, MSc., RMCCM

Microbiology Lab Manager

EMSL Canada Inc., Mississauga Laboratory



# FINAL REPORT

Assessment of Antimicrobial Activity

Using a Time-Kill Procedure

Order Number: 551801924

## PREPARED FOR


Franke Kindred Canada | 1000 Franke Kindred Road  
Midland, Ontario L4R 4K9

Sneha Panchal, M.Sc., RMCCM

Laboratory Director

2/28/2018

EMSL Canada Inc.  
2756 Slough Street, Mississauga, ON L4T 1G3  
Phone: (289) 997-4602 Fax: (289) 997-4607 Web: [www.emsl.com](http://www.emsl.com)

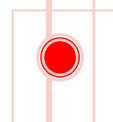

## CERTIFICATE OF ANALYSIS

---

**CLIENT:** FRANKE KINDRED CANADA

**PRODUCT:** OZONATED TAP WATER

**CONTACT:** CORY MACEY

**SAMPLE RECEIVED:** 02/01/2018

**PROJECT:**

**REPORT DATE:** 02/28/2018

ASSESSMENT OF ANTIMICROBIAL  
ACTIVITY USING ASTM METHOD  
E2315

**CHALLENGE BACTERIA:**

*Clostridium difficile* spores

*Legionella pneumophila*

### I. EXPERIMENTAL SUMMARY

The testing procedure was designed after discussions between EMSL Canada Inc. and Franke Kindred Canada. The procedure is based on ASTM E2315 method guidelines and conducted on ozonated water samples of concentrations 2.0 ppm to demonstrate its effectiveness at killing *Clostridium difficile* (ATCC 9689) spores and *Legionella pneumophila*. The testing was conducted in the Mississauga Microbiology Laboratory.

### II. PROCEDURE

The testing was done to determine the effectiveness of ozonated water sample of concentration 2.0 ppm (provided by Cory Macey) at killing *Clostridium difficile* (ATCC 9689) spores for 20 seconds exposure time and *Legionella pneumophila* for 5 minutes exposure time.

**Culture preparation:**

*Clostridium difficile* spores: A pure culture of *Clostridium difficile* (ATCC 9689) was streaked on to Healthlink Anaerobic Blood Agar plates and incubated at 35°C for up to 10 days under anaerobic conditions. Following incubation, the surface of agar plate was scraped and spores were harvested. The control spore suspension was adjusted to a concentration of >90%.

*Legionella pneumophila*: A pure culture of *Legionella pneumophila* was streaked onto Buffered Charcoal Yeast Extract agar for 48 hours at 35°C. The control suspension of the test microorganism was standardized to a minimum concentration of  $1.0 \times 10^6$  CFU/mL.



Test and control substances were dispensed in identical volumes to sterile test tubes. Independently, Test and Control substances were inoculated with the test microorganism, mixed and incubated.

Control suspensions were immediately plated to represent the concentration present at the start of the test, or time zero.

At the conclusion of each contact time, a volume of the liquid test solution was neutralized. Dilutions of the neutralized test solution were plated on to appropriate agar plates and incubation temperatures to determine the surviving microorganisms at the respective contact times.

Reductions of microorganisms were calculated by comparing initial microbial concentrations to surviving microbial concentrations.

All tests were performed in duplicates and counts averaged.

#### Calculations:

Calculations were based on the following:

$$\text{Percentage reduction} = (B-A/B) \times 100$$

$$\text{Log 10 Reduction} = \text{Log (B/A)}$$

Where:

B = Number of viable test microorganisms in the control substance immediately after inoculation

A = Number of viable test microorganisms in the test substance after the contact time

### III. EXPERIMENTAL RESULTS

Test Microorganism	Test Substance	Contact Time	CFUs/mL	Percent Reduction Compared to Control at Time Zero	Log10Reduction compared to Control at Time Zero
<i>Legionella pneumophila</i>	CONTROL	Time Zero	2.0E+06	n/a	
	<b>SAMPLE A:</b> 2.0 PPM	5 minutes	<1.00E+01	99.99%	5.3



Test Microorganism	Test Substance	Contact Time	CFUs/mL	Percent Reduction Compared to Control at Time Zero	Log10 Reduction compared to Control at Time Zero
<i>C.difficile</i> spores	CONTROL	Time Zero	3.9E+03	n/a	
	<b>SAMPLE A:</b> 2.0 PPM CONC.	20 seconds	1.4E+03	64.10%	0.44

**IV. CONCLUSIONS/OBSERVATIONS**

Sample A (concentration 2.0 ppm) caused a 99.99% reduction on *L.pneumophila* upon 5 minutes exposure.



Sneha Panchal, MSc., RMCCM

Laboratory Director

EMSL Canada Inc., Mississauga Laboratory

